
growth curve batch 3 at 35C

Growth curve of 32 strains

| | genome_ie | d exp_id | genome_ | _id expid | genom | e_id exp_id | | genon | ne_id exp_id |
|---|-----------|-----------|---------|-----------|-------|-------------|----|-------|--------------|
| 1 | g2 | H2M3R1 9 | g11 | L2M4R1 17 | g22 | fp1-2 | 25 | g30 | gp1-3 |
| 2 | g3 | H2M3R2 10 | g13 | L3M5R1 18 | g23 | fp1-3 | 26 | g31 | bg-2 |
| 3 | g4 | H3M1R1 11 | g15 | L4M2R2 19 | g24 | fp2-2 | 27 | g32 | bg-3 |
| 4 | g5 | H3M3R2 12 | g16 | L4M3R3 20 | g25 | crp1-2 | 28 | g33 | pms-1 |
| 5 | g6 | H3M4R1 13 | g17 | L4M4R1 21 | g26 | crp1-3 | 29 | g34 | pms-2 |
| 6 | g8 | H4M5R1 14 | g19 | L4M7R1 22 | g27 | crp2-2 | 30 | g35 | pms-3 |
| 7 | g9 | L1M2R2 15 | g20 | L3M1R1 23 | g28 | gp1-1 | 31 | g36 | ppf-1 |
| 8 | g10 | L2M2R1 16 | g21 | src-2 24 | g29 | gp1-2 | 32 | g37 | 40 th- 1 |

- revive glycerol stock of rhizobia on TY agar
 - prepare labelled plates
 - use a sterilized inoculation loop to streak the glycerol stock on TY agar
 - incubate the petri plates at 30C for 2-3 days until the colony is identifiably large
- prepare ensifer inoculum
 - prepare a DW96 filled with 300uL of PBS
 - pick colonies from the petri plate and dissolve the colonies in the PBS. Use wells in column 1, 5, and 9
 - dispense 100 uL of the dissolved colony to a NUNC96 and measure od at 600 nm
 - standardize the OD600 to 0.1. For target od = 0.1 od_target and target inoculum volume = 100 uL volume_target, and given the measured od of dissolved colonies od_dissolved, use the code below to calculate the needed volume of the dissolved colonies volume_dissolved and the pbs volume volume pbs
 - prepare another NUNC96, follow the table to add required volume of dissolved colonies and pbs to make 100 uL of standardized inoculum
 - measure the od of the standardized inoculum as od_inoculum
- prepare 100 uL of TY medium in NUNC96
 - inoculate 2 uL of the standardized inoculum to each well

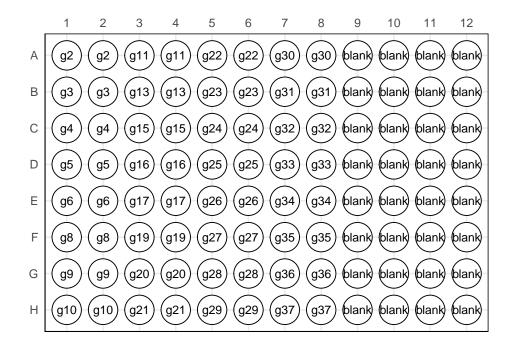

| genome_id | well | od_dissolved | od_target | volume_target | volume_dissolved | volume_pbs |
|-----------|------|--------------|-----------|---------------|------------------|------------|
| g2 | A01 | 0.449125 | 0.1 | 100 | 22.27 | 77.73 |
| g3 | B01 | 0.844125 | 0.1 | 100 | 11.85 | 88.15 |
| g4 | C01 | 0.340125 | 0.1 | 100 | 29.40 | 70.60 |
| g5 | D01 | 0.221125 | 0.1 | 100 | 45.22 | 54.78 |
| g6 | E01 | 0.227125 | 0.1 | 100 | 44.03 | 55.97 |
| g8 | F01 | 0.330125 | 0.1 | 100 | 30.29 | 69.71 |
| g9 | G01 | 0.214125 | 0.1 | 100 | 46.70 | 53.30 |
| g10 | H01 | 0.398125 | 0.1 | 100 | 25.12 | 74.88 |
| g11 | A03 | 0.245125 | 0.1 | 100 | 40.80 | 59.20 |
| g13 | B03 | 0.115125 | 0.1 | 100 | 86.86 | 13.14 |
| g15 | C03 | 0.817125 | 0.1 | 100 | 12.24 | 87.76 |
| g16 | D03 | 0.190125 | 0.1 | 100 | 52.60 | 47.40 |
| g17 | E03 | 0.209125 | 0.1 | 100 | 47.82 | 52.18 |
| g19 | F03 | 0.104125 | 0.1 | 100 | 96.04 | 3.96 |
| g20 | G03 | 0.169125 | 0.1 | 100 | 59.13 | 40.87 |
| g21 | H03 | 0.196125 | 0.1 | 100 | 50.99 | 49.01 |
| g22 | A05 | 0.474125 | 0.1 | 100 | 21.09 | 78.91 |
| g23 | B05 | 0.384125 | 0.1 | 100 | 26.03 | 73.97 |
| g24 | C05 | 0.302125 | 0.1 | 100 | 33.10 | 66.90 |
| g25 | D05 | 0.141125 | 0.1 | 100 | 70.86 | 29.14 |
| g26 | E05 | 0.304125 | 0.1 | 100 | 32.88 | 67.12 |
| g27 | F05 | 0.319125 | 0.1 | 100 | 31.34 | 68.66 |
| g28 | G05 | 0.178125 | 0.1 | 100 | 56.14 | 43.86 |
| g29 | H05 | 0.194125 | 0.1 | 100 | 51.51 | 48.49 |
| g30 | A07 | 0.318125 | 0.1 | 100 | 31.43 | 68.57 |
| g31 | B07 | 0.255125 | 0.1 | 100 | 39.20 | 60.80 |
| g32 | C07 | 0.557125 | 0.1 | 100 | 17.95 | 82.05 |
| g33 | D07 | 0.317125 | 0.1 | 100 | 31.53 | 68.47 |
| g34 | E07 | 0.292125 | 0.1 | 100 | 34.23 | 65.77 |
| g35 | F07 | 0.336125 | 0.1 | 100 | 29.75 | 70.25 |
| g36 | G07 | 0.174125 | 0.1 | 100 | 57.43 | 42.57 |
| g37 | H07 | 0.403125 | 0.1 | 100 | 24.81 | 75.19 |