
growth curve batch 5, a second run at 25C

Growth curve of 32 strains

	genome_id	exp_id	genome_id	exp_id		genome_id	exp_id		genome_id	exp_id
1	g2	H2M3R1 9	g11	L2M4R1	17	g22	fp1-2	25	g30	gp1-3
2	g3	H2M3R2 10	g13	L3M5R1	18	g23	fp1-3	26	g31	bg-2
3	g4	H3M1R1 11	g15	L4M2R2	19	g24	fp2-2	27	g32	bg-3
4	g5	H3M3R2 12	g16	L4M3R3	20	g25	crp1-2	28	g33	pms-1
5	g6	H3M4R1 13	g17	L4M4R1	21	g26	crp1-3	29	g34	pms-2
6	g8	H4M5R1 14	g19	L4M7R1	22	g27	crp2-2	30	g35	pms-3
7	g9	L1M2R2 15	g20	L3M1R1	23	g28	gp1-1	31	g36	ppf-1
8	g10	L2M2R1 16	g21	$\operatorname{src-2}$	24	g29	gp1-2	32	g37	40 th- 1

- revive glycerol stock of rhizobia on TY agar
 - prepare labelled plates
 - use a sterilized inoculation loop to streak the glycerol stock on TY agar
 - incubate the petri plates at 30C for 2-3 days until the colony is identifiably large
- prepare ensifer inoculum
 - prepare a DW96 filled with 300uL of PBS
 - pick colonies from the petri plate and dissolve the colonies in the PBS. Use wells in column 1, 5, and 9
 - dispense 100 uL of the dissolved colony to a NUNC96 and measure od at 600 nm
 - standardize the OD600 to 0.1. For target od = 0.1 od_target and target inoculum volume = 100 uL volume_target, and given the measured od of dissolved colonies od_dissolved, use the code below to calculate the needed volume of the dissolved colonies volume_dissolved and the pbs volume volume pbs
 - prepare another NUNC96, follow the table to add required volume of dissolved colonies and pbs to make 100 uL of standardized inoculum
 - measure the od of the standardized inoculum as od_inoculum
- prepare 100 uL of TY medium in NUNC96
 - inoculate 2 uL of the standardized inoculum to each well

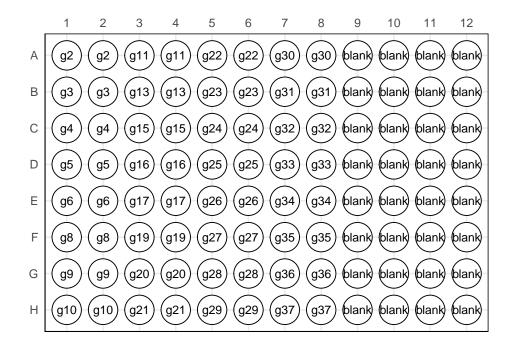

genome_id	well	od_dissolved	od_target	volume_target	volume_dissolved	volume_pbs
g2	A01	0.511125	0.1	100	19.56	80.44
g3	B01	0.572125	0.1	100	17.48	82.52
g4	C01	0.537125	0.1	100	18.62	81.38
g5	D01	0.531125	0.1	100	18.83	81.17
g6	E01	0.260125	0.1	100	38.44	61.56
g8	F01	0.703125	0.1	100	14.22	85.78
g9	G01	0.416125	0.1	100	24.03	75.97
g10	H01	0.418125	0.1	100	23.92	76.08
g11	A03	0.210125	0.1	100	47.59	52.41
g13	B03	0.275125	0.1	100	36.35	63.65
g15	C03	0.712125	0.1	100	14.04	85.96
g16	D03	0.548125	0.1	100	18.24	81.76
g17	E03	0.335125	0.1	100	29.84	70.16
g19	F03	0.272125	0.1	100	36.75	63.25
g20	G03	0.294125	0.1	100	34.00	66.00
g21	H03	0.315125	0.1	100	31.73	68.27
g22	A05	0.313125	0.1	100	31.94	68.06
g23	B05	0.291125	0.1	100	34.35	65.65
g24	C05	0.560125	0.1	100	17.85	82.15
g25	D05	0.451125	0.1	100	22.17	77.83
g26	E05	0.610125	0.1	100	16.39	83.61
g27	F05	0.310125	0.1	100	32.25	67.75
g28	G05	0.285125	0.1	100	35.07	64.93
g29	H05	0.283125	0.1	100	35.32	64.68
g30	A07	0.615125	0.1	100	16.26	83.74
g31	B07	0.371125	0.1	100	26.95	73.05
g32	C07	0.315125	0.1	100	31.73	68.27
g33	D07	0.227125	0.1	100	44.03	55.97
g34	E07	0.439125	0.1	100	22.77	77.23
g35	F07	0.328125	0.1	100	30.48	69.52
g36	G07	0.344125	0.1	100	29.06	70.94
g37	H07	0.322125	0.1	100	31.04	68.96